



Custom Form Designer Store and Forward Telemedicine Software

PROFESSIONAL


Second Opinion Professional allows you to capture high resolution color or grayscale images, video and sound clips, and other medical information and store them in a database folder for review. Using both real time and store and forward technologies, you can send the information to remote locations for a second opinion.

The Custom Form Designer program allows users to create their own custom forms by simply dragging and dropping fields on to a page. There are virtually no limits to the number of pages or the number of fields you can have. Once finished with the design process, all you have to do is generate the form, decide the subject and type of form and the industry standard xBase databases, and screen definitions are automatically created. These are not word processing or email forms; they are true database entry applications. Data that is collected can be reported on with just about any report generator program on the market.


The Form Properties can be set to read-write or read-only, and free or restricted distribution. Groups can create custom forms to be fully shared, and be assured that the design of the form does not get altered.

When a Custom Form is transmitted to another site, the form definition travels with the data. That means the person receiving the data does not have to have the form ahead of time. It will be automatically added to their system.


● Design

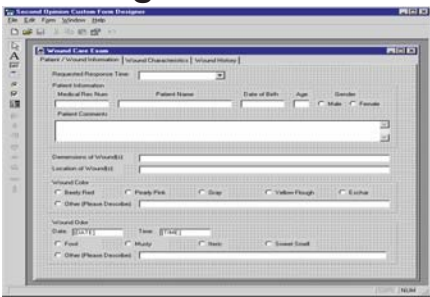

● Generate

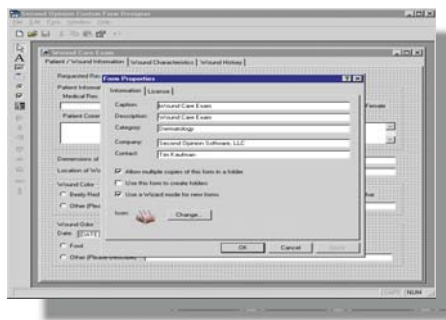

● Execute

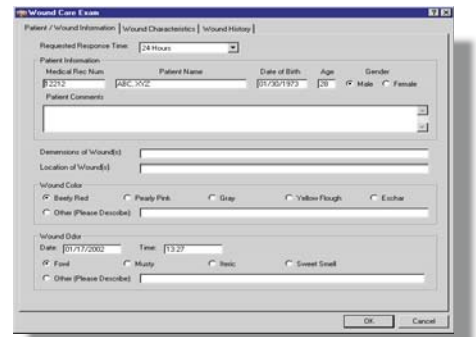

Second Opinion Telemedicine Solutions, Inc.
3848 Del Amo Blvd, Suite 303
Torrance, CA 90503
(310) 802-6300
www.soTelemed.com
soSales@soTelemed.com