



Store and Forward Telemedicine Software

PROFESSIONAL

This telehealth application allows you to capture high resolution color or grayscale images, video and sound clips, and other medical information and store them in a database folder for review. Using both real time and store and forward technologies, you can send the information to remote locations for a second opinion.


Second Opinion™ Professional can be used in all medical specialties including Cardiology, Ophthalmology, Dermatology, Pathology, Gynecology and Forensic Science applications.

Features:


- Standalone Client Software
- Web Client / Server Software
- Network Client / Server Software
- Compatible w/ All On-Line Image Libraries
- Compatible w/ All On-Line Image Storage Systems

Free Receive / Viewer Application Available

• Organize

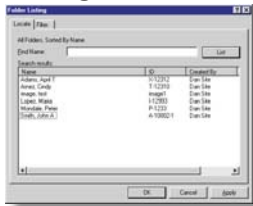

• Manage Your Documents

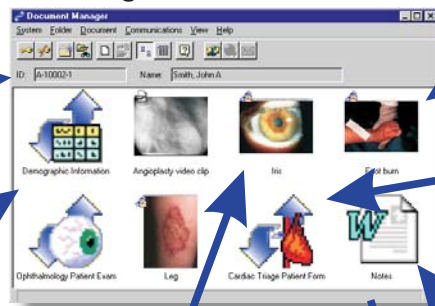


• Review



• Follow Protocols




• Gather Information

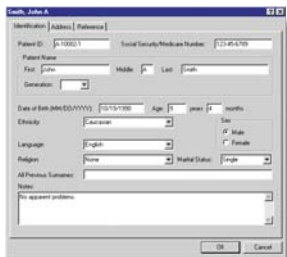

• Capture Pictures



• Share Data

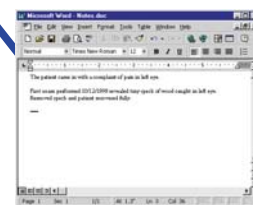

• Link to Other Sites




Second Opinion Dramatically Improves the Way Information is Shared Within the Medical Community


Second Opinion Professional

The **Second Opinion™ Form Designer** program allows users to create their own custom forms by simply placing fields on to a page. There are virtually no limits to the number of pages or the number of fields you can have. Once finished with the design process, all you have to do is generate the form, decide the subject and type of form and the industry standard xBase databases, and screen definitions are automatically created. These are not word processing or email forms; they are true database entry applications. Data that is collected can be reported on with the **Second Opinion™ Report Designer** or just about any report generator program on the market.



• Design

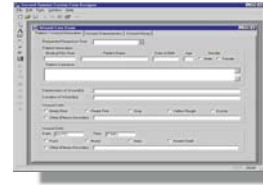

• Execute

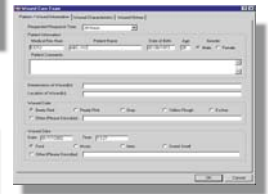

• Generate

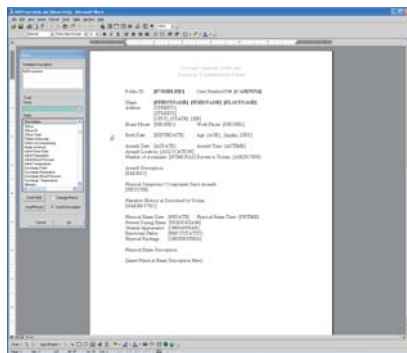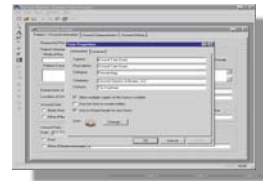

The **Second Opinion™ Report Designer** transforms Microsoft Word into a report writer. Templates are created so data from within **Second Opinion™ Professional** can be transferred to the Word templates. Once this is done, the content can be, adjusted and new narrative added. Pictures from within the patients folder can be added as well and printed as a final report. There is no limit to the size of a report.

The **Second Opinion™ Professional Software Developer Kit (SDK)** allows developers to create custom add-ins. It is also useful for developers who wish to add video capture, scanning, imaging, communications, user security, and much more to their existing applications. **Second Opinion™ Professional** acts as an OLE Automation Server and can be run as a front-end application fully visible or can be executed hidden as a back-end process.

• Sample Code



• Objects



• Methods & Properties



All of the Methods, Properties, and Sequence of Events throughout the program set are fully exposed and documented, making the creation and integration tasks very straight forward.



Second Opinion

Telemedicine Solutions, Inc.

Second Opinion Telemedicine Solutions, Inc.

3848 Del Amo Blvd, Suite 303

Torrance, CA 90503

(310) 802-6300

www.soTelemed.com

soSales@soTelemed.com